



High-Throughput SPE/GC-MS Method for Quantification of Methamphetamine and Amphetamine In Urine Samples of Drug Users

Sumate Thiangthum, Orapin Tanunkat

Bureau of Drug and Narcotic, Department of Medical Sciences, Ministry of Public Health, Nonthaburi 11000, Thailand.

Introduction

Methamphetamine is the most common drug of abuse in Thailand. Bureau of Drug and Narcotic, Department of Medical Sciences, as the national drug testing laboratory, has analyzed more than hundred thousand urine specimens collected from drug users each year. It is necessary to develop a high-throughput method capable of simultaneously determining methamphetamine and its metabolite, amphetamine, in urine specimens.

Method

A simple and rapid GC-MS method was developed using automated solid-phase extraction for sample clean up. Phentermine was used as internal standard. Derivatization was performed with pentafluoropropionic acid at 65°C for 25 minutes. Chromatography was conducted on a fused silica capillary column and analytes were determined in selected-ion-monitoring (SIM) mode. Mass spectra of pentafluoropropionyl derivatives showed peak at m/z 190, 118 and 91 for amphetamine, at m/z 204, 160 and 118 for methamphetamine and at m/z 204, 132 and 91 for phentermine. The method was fully validated according to the current recommendations of the USFDA bioanalytical method validation guidance.



Table 1 Precision, accuracy and recovery data

Analyte	Concentration (ng/mL)	Intraday		Inter-day		%Recovery
		%CV	%RD	%CV	%RD	
Methamphetamine	100	0.87	8.53	12.98	3.78	-
	250	3.92	4.04	9.12	3.37	91.80
	2,500	1.54	2.39	5.68	-0.56	76.89
	4,000	6.30	8.54	3.08	0.13	90.86
Amphetamine	100	4.15	16.76	3.84	17.23	-
	250	10.17	-0.54	9.52	3.52	99.46
	2,500	8.42	5.98	10.04	-0.34	75.58
	4,000	6.14	8.33	5.23	9.08	75.40

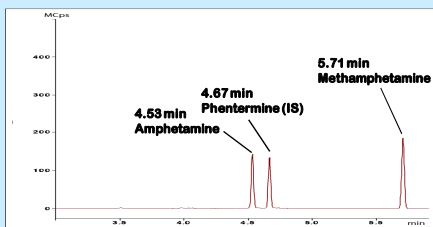


Fig. 1 Chromatogram of a spiked standard mixture containing methamphetamine (2,500 ng/mL), amphetamine (2,500 ng/mL), and phentermine (1,000 ng/mL)

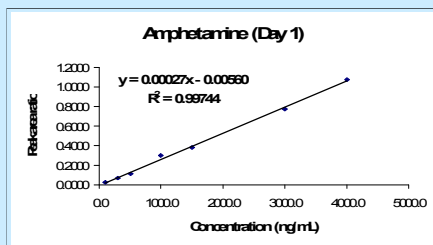


Fig. 2 Calibration curve of amphetamine spiked standard in range of 100 -4,000 ng/mL

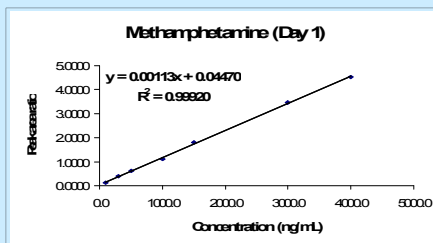


Fig. 3 Calibration curve of methamphetamine spiked standard in range of 100 -4,000 ng/mL

Results

There were no interfering peaks from endogenous components in blank urine chromatograms. Calibration curves were linear cover 100-4,000 ng/mL with correlation coefficients greater than 0.992 (Fig.2-3). Recovery was 75.4-99.5% for amphetamine and 76.9-91.8% for methamphetamine (Table 1). Accuracy and precision were performed at 4 different concentrations cover the calibration range. Accuracy, expressed as relative deviation (%RD) was less than 16.8% for intra-day and less than 17.2% for inter-day (Table 1). Precision, expressed as coefficient of variation (%CV) was less than 10.2% for intra-day and less than 13.0% for inter-day (Table 1). Stability of amphetamine and methamphetamine in various storage conditions was also determined.

Conclusion

The method was proved to be high-throughput and successfully used in our laboratory to quantify amphetamine and methamphetamine in urine samples of drug users.

Acknowledgment

The financial support from Department of Medical Sciences foundation was greatly appreciated.

Literature cited

- U.S. Department of Health and Human Services, Food and Drug Administration, Center for Drug Evaluation and Research (CDER), Center for Veterinary Medicine (CVM). Guidance for Industry: Bioanalytical Method Validation. Internet at <http://www.fda.gov/cder/guidance/index.htm> (2001).
- Whitmire M, Ammerman J, de Liso P, Killmer J, Kyle D, Mainstone E, Porter L, Zhang T. LC-MS/MS bioanalysis method development. Validation, and sample analysis: Points to consider when conducting nonclinical and clinical studies in accordance with current regulatory guidances. J Anal. Bioanal Techniques 2011;54:1-10. Internet at <http://dx.doi.org/10.4172/2155-9872.54-001>.
- Goldberger BA, Huestis MA, Wilkins DG. Commonly practiced quality control and quality assurance procedures for gas chromatography/mass spectrometry analysis in forensic urine drug-testing laboratories. Forensic Sci Rev 1997;9(2):60-80.
- Viswanathan CT, Bansal S, Booth B, DeStefano AJ, Rose MJ, Sallstad J, Shah VP, Skelly JP, Swann PG, Weiner. Workshop/conference report-Quantitative bioanalytical methods validation and implementation: best practices for chromatographic and ligand binding assays. The AASP Journal 2007;9(1):E30-42.
- Bansal S, DeStefano A. Key elements of bioanalytical method validation for small molecules. The AASP Journal 2007;9(1):E109-14.